


CAPP *Rondo*

TUBE ROTATOR

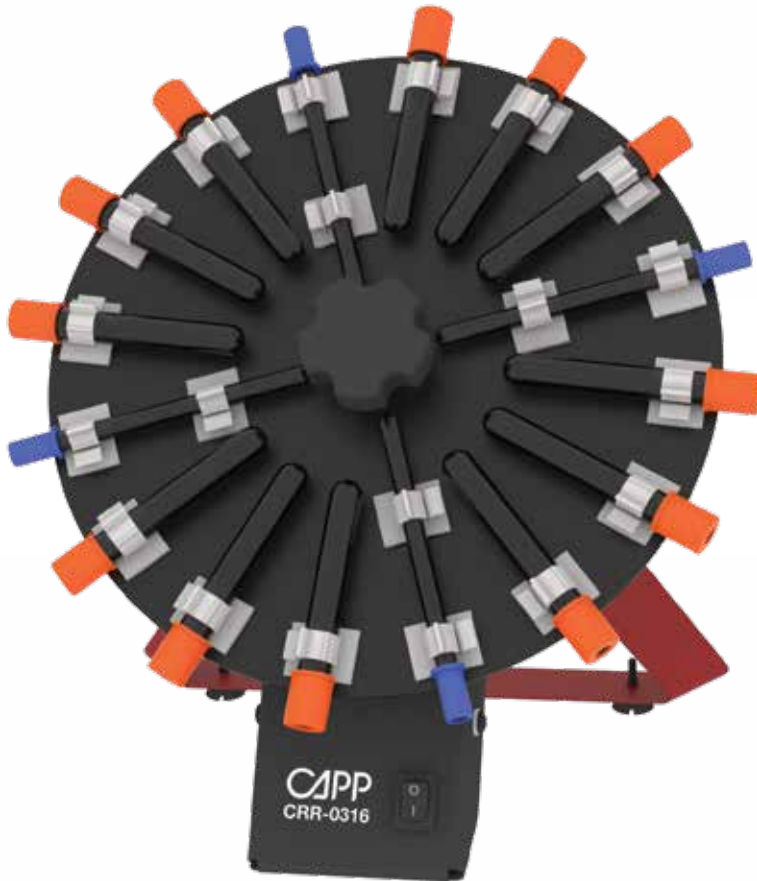

CAPP *Rondo* Tube Rotator has a strong, metal design which, along with the rubber feet, ensures its total stability during operation.

The compact design uses a minimal amount of benchtop space.

The gentle but effective mixing occurs at a fixed speed of 30 RPM with an accuracy of +/- 2 RPM.

CAPP *Rondo* Tube Rotator has a maximum capacity of 16 tubes, consisting of 4 ESR and 12 blood collection tubes. Its maximum operational weight capacity is 1.5kg. The disc plate is fixed at an angle of 38 degrees and is easily removed for cleaning or replacement.

CAPP *Rondo* Tube Rotator operates silently in the continuous mode.

Specifications:

Max Capacity	16 Tubes (4 ESR Tubes and 12 Blood Collection Tubes)
Motor Type	AC Synchronous
Speed	Fixed Speed of 30 RPM
Operation Mode	Continuous / Infinite
Speed Accuracy	+/- 2 RPM
Angle Of Disc Plate	Fixed angle of 38 degrees
Noise Level	Noiseless operation
Load Capacity	1.5kg
Dimensions	278x183x252mm
Protection Class	IP 21

Ordering Information

Cat No.	Description
CRR-0316	CAPPRondo Tube Rotator for 16 blood tubes
CRR-0316-110	CAPPRondo Tube Rotator for 16 blood tubes, 110V

