

EXPell Cryo Tubes


Superior by Design™



CAPP offers a unique portfolio of cryo tubes in 5 sizes: 0.5mL, 1.0mL, 2.0mL, 4.0mL and 5.0mL.

EXPell Cryo Tubes are designed with a focus on the best user experience and convenience. Featured with single-handed QuickSeal™ closure, the external thread design allows tightening or removal with just half a turn. These are one of the only cryo tubes with external thread and star foot which can be stored in 10x10 cryoboxes, saving up to 23% storage capacity compared to conventional cryo solutions. They are also compatible with common workstations on the market. Both the caps and the cryo tubes are made of specially formulated polypropylene, having the same coefficient of expansion and therefore enhancing the leakproof qualities at various temperatures.

EXPell Cryo Tubes come packed in resealable stand-up pouches, saving benchtop space. They are suitable for storage in the vapour phase of liquid nitrogen.

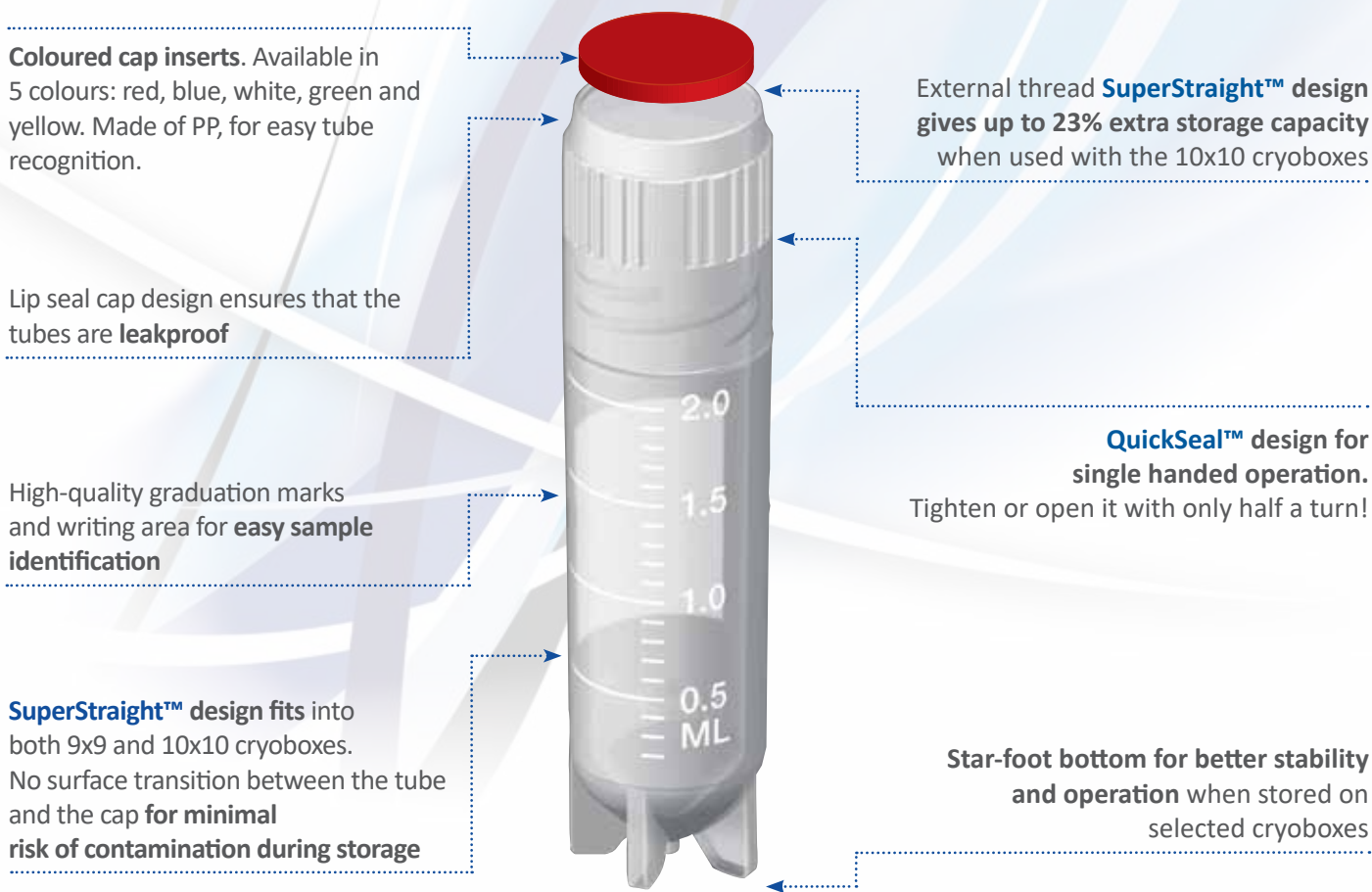

Ordering Information

EXPell Cryo Tubes

Cat. No.	Description
5040502C	Expell Cryo Tube 0,5mL, sterile, w/o print, bulk, case with 2x500 pcs (tubes) and 1x1000 pcs (caps)
5041000C	Expell Cryo Tube 1.0mL, 10x100 pcs.
5042000C	Expell Cryo Tube 2.0mL, 10x100 pcs.
5042001C	Expell Cryo Tube 2.0mL, non-sterile, bulk, case with 2x500 pcs (tubes) and 1x1000 pcs (caps)
5044000C	Expell Cryo Tube 4.0mL, 5x100 pcs.
5045000C	Expell Cryo Tube 5.0mL, 5x100 pcs.

Coloured Cap Inserts for EXPell Cryo Tubes

Cat. No.	Description
504R000C	Red Cap Insert for Expell Cryo Tube, 1x500 pcs.
504B000C	Blue Cap Insert for Expell Cryo Tube, 1x500 pcs.
504G000C	Green Cap Insert for Expell Cryo Tube, 1x500 pcs.
504W000C	White Cap Insert for Expell Cryo Tube, 1x500 pcs.
504Y000C	Yellow Cap Insert for Expell Cryo Tube, 1x500 pcs.
MLD-50	Cryo Tubes Decapping Tool, 50 pcs.